



Features:

- Flow rate exceeds the requirements of NFPA and CSA
- Accepts quick disconnect gas specific adapters, or DISS adapters depending on which front assembly is ordered
- 100 % Hydrostatic testing of each outlet and cleaned for medical gas service
- Indexed to eliminate interchangeability of gas services
- Up to 1" [25mm] wall thickness adjustment
- 360° Swivel Inlet Pipe for easy installation
- Universal Rough-In Assembly to accept Quick Disconnects or DISS front adapters (Can be interchanged at any time)
- Outlets are also available for M.R.I. use

General Specifications:

The medical gas outlet shall be an Amico Alert-1 series.

Outlets shall be manufactured with a 6-1/2" [165 mm] length type "K" 1/2" [12.7 mm] outside diameter (3/8" nominal) size copper inlet pipe stub, which is silver brazed to the outlet body. The body shall be of 1-5/16" [33 mm] diameter, one piece brass construction. For positive pressure gas services, the outlet shall be equipped with a primary and secondary check valve and the secondary check valve shall be rated at a maximum of 200 psi [1,379 kPa] in the event the primary check valve is removed for maintenance. Outlet bodies shall be gas specific by indexing each gas service to a gas specific dual pin indexing arrangement on the respective identification module.

A large color coded front plate shall be used for ease of gas identification and aesthetic appeal.

With the back rough-in mounted the outlet shall adjust from 3/8" [10mm] up to 1" [25mm] variation in wall thickness.

The latch-valve assembly shall be Ohmeda Diamond™ Quick Disconnect, Chemetron Quick Disconnect, Puritan-Bennett Quick Disconnect, or D.I.S.S. (Diameter Index Safety System) depending on the type ordered.

All M.R.I. outlets are manufactured from nonferrous materials.

All outlets shall be cleaned and degreased for medical gas service, factory assembled and tested.


The outlets are UL Listed and ETL Listed.

Amico products comply with NFPA-99 and CSA.

Technical Specifications

NOTE:
Amico recommends that you have a blank space between outlets and slides.
Maximum 3 outlets may be ganged together without additional support.

Inch
[mm]

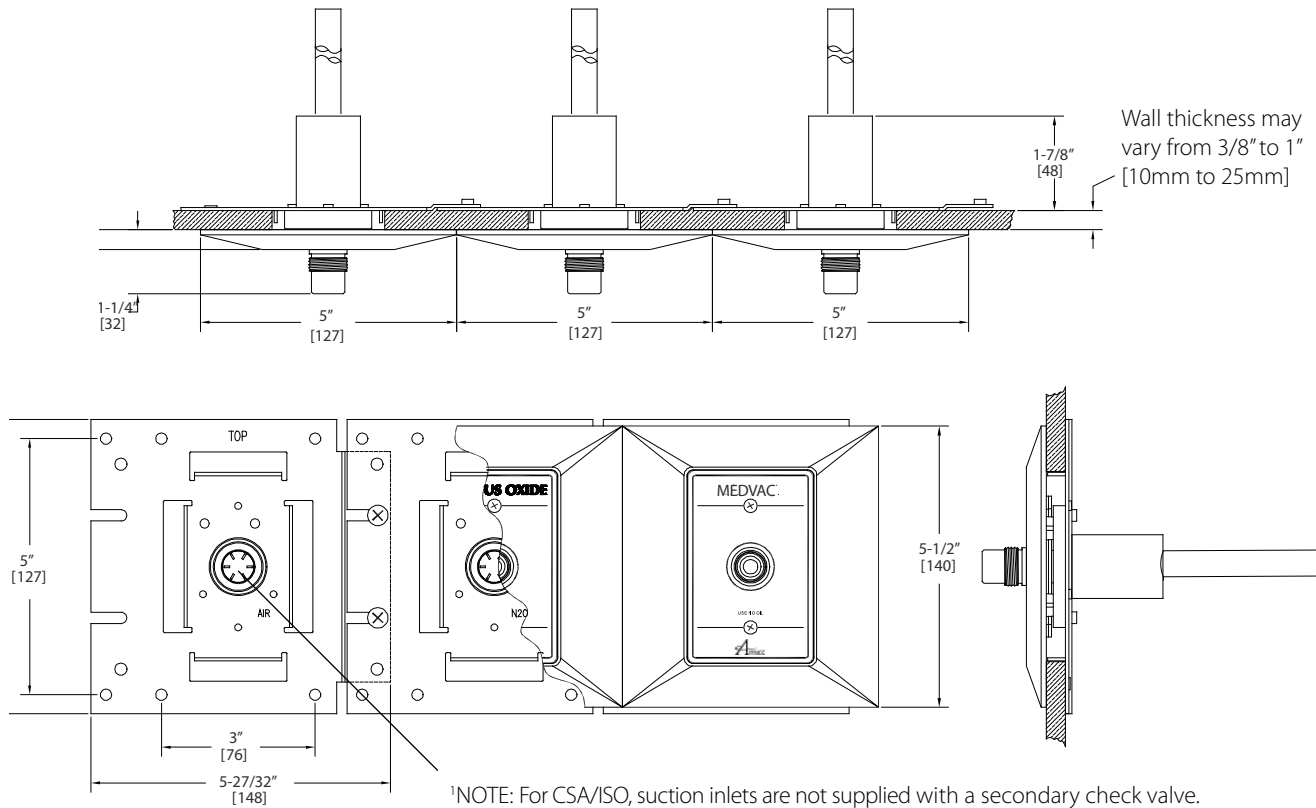

¹NOTE: For CSA/ISO, suction inlets are not supplied with a secondary check valve.
A pressure plug is provided for testing purposes, rated at 150 psi [1,034 kPa].

Model Numbers

O-DISCEI-L-XXX

DIS = DISS
QD = Ohmeda Compatible
CH = Chemetron Compatible
PB = Puritan-Bennett

The "L" defines the Language:

U = English (NFPA) Oxy, Air, Vac and WAG only
E = English (CSA/ISO)
F = French (CSA/ISO)
S = Spanish (NFPA)

The "XXX" defines the Gas:

Oxygen	=	OXY
Medical Air	=	AIR
Nitrous Oxide	=	N2O
*Nitrogen	=	NIT
**Carbon Dioxide	=	CO2
MedVac ¹	=	VAC
Waste Anesthetic Gas Disposal (NFPA)	=	WAG
Anesthetic Gas Scavenging System ¹ (CSA)	=	AGS
Instrument Air	=	INA

* = Available in DISS only

** = Not available for Puritan-Bennett

¹ Please see note above

Example: Oxygen NFPA - DISS, Oxygen = **O-DISCEI-U-OXY**

Amico Corporation 85 Fulton Way, Richmond Hill, ON L4B 2N4, Canada | 71 East Industry Court, Deer Park, NY 11729, USA
Toll Free Tel: 1.877.462.6426 | Toll Free Fax: 1.866.440.4986 Tel: 905.764.0800 | Fax: 905.764.0862