


## Application

- Supports the whole body including the head
- Can be used to lift patients who have little control over their head and body
- Can be used to lift and transfer from all patient surfaces, including the floor

**Maximum load = 1000 lbs**

Part Number	Material	Size	Head Support
SLG-GO-BASIC-WV-HS-S	Polyester	Small (S)	Yes
SLG-GO-BASIC-WV-HS-M	Polyester	Medium (M)	Yes
SLG-GO-BASIC-WV-HS-L	Polyester	Large (L)	Yes
SLG-GO-BASIC-WV-HS-XL	Polyester	Extra Large (XL)	Yes
SLG-GO-BASIC-MSH-HS-S	Mesh	Small (S)	Yes
SLG-GO-BASIC-MSH-HS-M	Mesh	Medium (M)	Yes
SLG-GO-BASIC-MSH-HS-L	Mesh	Large (L)	Yes
SLG-GO-BASIC-MSH-HS-XL	Mesh	Extra Large (XL)	Yes
SLG-GO-BASIC-WV-S	Polyester	Small (S)	No
SLG-GO-BASIC-WV-M	Polyester	Medium (M)	No
SLG-GO-BASIC-WV-L	Polyester	Large (L)	No
SLG-GO-BASIC-WV-XL	Polyester	Extra Large (XL)	No
SLG-GO-BASIC-MSH-S	Mesh	Small (S)	No
SLG-GO-BASIC-MSH-M	Mesh	Medium (M)	No
SLG-GO-BASIC-MSH-L	Mesh	Large (L)	No
SLG-GO-BASIC-MSH-XL	Mesh	Extra Large (XL)	No

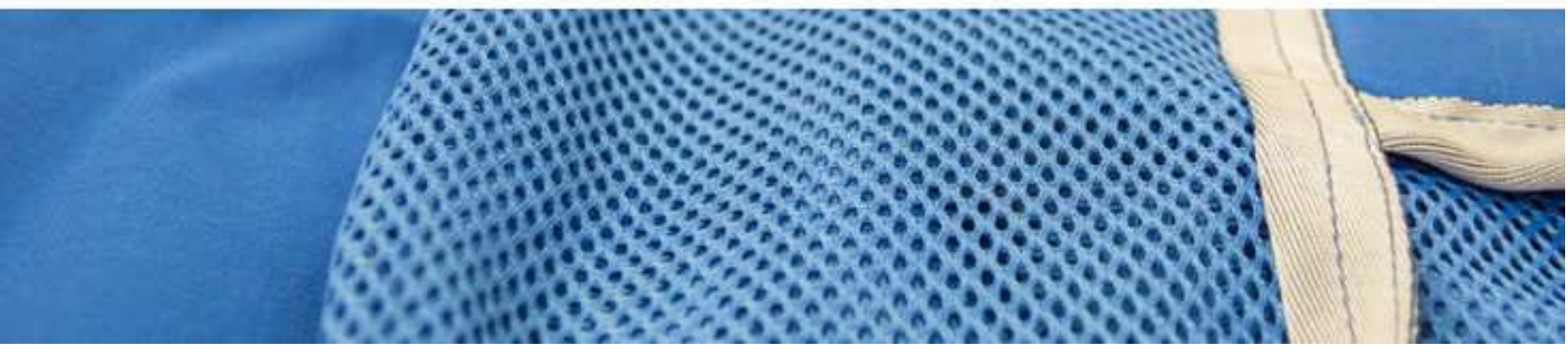

# GoLong

## Application

- Can be used to reposition the patient in bed
- Ideal for transferring patients to other surfaces including beds and stretchers
- Can be left underneath the patient

Maximum load = 1000 lbs

Part Number	Material	Size
SLG-GO-LONG-WV-42X78	Polyester	42" x 78"
SLG-GO-LONG-WV-56X78	Polyester	56" x 78"
SLG-GO-LONG-MSH-42X78	Mesh	42" x 78"
SLG-GO-LONG-MSH-56X78	Mesh	56" x 78"



# GoHygiene

## Application

- Suitable for patients who have head control but reduced upper body strength
- Support around the body can be adjusted with the belt
- Large opening at the bottom allows access for toileting and hygiene

Maximum load = 1000 lbs

Part Number	Material	Size
SLG-GO-HYGIENE-WV-S	Polyester	Small (S)
SLG-GO-HYGIENE-WV-M	Polyester	Medium (M)
SLG-GO-HYGIENE-WV-L	Polyester	Large (L)
SLG-GO-HYGIENE-WV-XL	Polyester	Extra Large (XL)



# GoComfort

## Application

- Supports the whole body except the head
- Can be used to provide good support for amputees
- Can be used to lift patients who have little control over their body
- Can lift from a seated or laying position
- Patients that require additional support around hip and thighs

Maximum load = 1000 lbs

Part Number	Material	Size
SLG-GO-COMFORT-WV-S	Polyester	Small (S)
SLG-GO-COMFORT-WV-M	Polyester	Medium (M)
SLG-GO-COMFORT-WV-L	Polyester	Large (L)
SLG-GO-COMFORT-WV-XL	Polyester	Extra Large (XL)



# GoCustom

All Slings are manufactured in-house and can be customized on demand. For your convenience, you may choose from a variety of options to get started...

Part Number	Type	Material	Size	Head Support
SLG-GO-CUSTOM	GoBasic	Polyester or Mesh	Small (S)	Yes or No
	GoLong	Polyester or Mesh	Medium (M)	Yes or No
	GoHygiene	Polyester or Mesh	Large (L)	Yes or No
	GoComfort	Polyester or Mesh	Extra Large (XL)	Yes or No



# Sizing Chart for Slings

	Slings	Material	S	M	L	XL
GoBasic	GoBasic With Head Support	Polyester	SLG-GO-BASIC-WV+HS-S	SLG-GO-BASIC-WV+HS-M	SLG-GO-BASIC-WV+HS-L	SLG-GO-BASIC-WV+HS-XL
	GoBasic With Head Support	Mesh	SLG-GO-BASIC-MSH+HS-S	SLG-GO-BASIC-MSH+HS-M	SLG-GO-BASIC-MSH+HS-L	SLG-GO-BASIC-MSH+HS-XL
	GoBasic No Head Support	Polyester	SLG-GO-BASIC-WV-S	SLG-GO-BASIC-WV-M	SLG-GO-BASIC-WV-L	SLG-GO-BASIC-WV-XL
	GoBasic No Head Support	Mesh	SLG-GO-BASIC-MSH-S	SLG-GO-BASIC-MSH-M	SLG-GO-BASIC-MSH-L	SLG-GO-BASIC-MSH-XL
GoLong	GoLong	Polyester	SLG-GO-LONG-WV-42X78		SLG-GO-LONG-WV-56X78	
	GoLong	Mesh	SLG-GO-LONG-MSH-42X78		SLG-GO-LONG-MSH-56X78	
GoHygiene	GoHygiene	Polyester	SLG-GO-HYGIENE-WV-S	SLG-GO-HYGIENE-WV-M	SLG-GO-HYGIENE-WV-L	SLG-GO-HYGIENE-WV-XL
	GoComfort	Polyester	SLG-GO-COMFORT-WV-S	SLG-GO-COMFORT-WV-M	SLG-GO-COMFORT-WV-L	SLG-GO-COMFORT-WV-XL

## Guidelines For Choosing Your Sling Size

	S	M	L	XL
Height	80-87 cm	85-92 cm	90-97 cm	95-102 cm
Width	36-40 cm	39-43 cm	42-45 cm	48-52 cm

## Washing Instructions



Normal washing at indicated temperature



Do not bleach



Tumble dry low



Do not iron



Do not dry clean

# Specifications & Listings

## Three GoLift Weight Capacities



Single Lift  
GO-LIFT-400      Lift Up to 400 lbs



Single Lift  
GO-LIFT-700      Lift Up to 700 lbs



Tandem Lifts<sup>1</sup>  
GO-LIFT-1000      Lift Up to 1000 lbs

## Safe Working Load (SWL)

400 lbs  
700 lbs  
1000 lbs<sup>1</sup>

## Dimensions

700/400 lbs:

- Length: 7 ¾"
- Width: 7 ¾"
- Height: 4 ¼"

1000 lbs:

- Length: 15 ½"
- Width: 7 ¾"
- Height: 4 ¼"
- Strap Length: 84"

## Unit Weight

400 lbs version: 8 lbs  
700 lbs version: 10 lbs  
1000 lb version: 20 lbs (Tandem Lifts)

## Lift Case

Flame retardant ABS

## Approvals

Certified to:  
Can/CSA-C22.2 No 60601-1:08, UL 60601-1:08  
Tested to: ISO 10535-06, CE.

## Safety

Emergency Stop  
Emergency Lowering Device  
Upper Limit Detection  
Lower Limit Detection  
Slack Tape Sensor  
Free Fall Brake (over speed governor)  
Low Battery and Dead Battery Alarms  
Soft Start and Stop  
Overload Protection  
Emergency Manual Lowering

## Hand Control

Capacitive touch  
Protection Class: IP67

## Service Life

10 years or 22,500 cycles

## Max Sound Level

<65 dB

## Max. Lifting Speed

No load: 2 inches/second  
150 lbs: 1.3 inches/second  
400 lbs: 1 inch/second  
700 lbs: 1 inch/second  
1000 lbs: 1 inch/second

## Batteries

High Capacity, Nickel Metal Hydride (Ni-Mh)  
Standard: 2 x 14.4V (2Ah)  
Optional: 2 x 14.4V (3.3 Ah) or (5 Ah)

## Battery Charger:

100-240 VAC 40 Watt

## Number of Lifts per Charge (Duty: 10/90)

375 with 185 lbs (25% of strap at midrange)<sup>2</sup>  
170 with 400 lbs (25% of strap at midrange)<sup>2</sup>  
102 with 700 lbs (25% of strap at midrange)<sup>2</sup>  
Charging time: 2-4 hours.

<sup>1</sup>Achieved using two GoLift 700 lifts in tandem

<sup>2</sup>Calculated using 5Ah Packs

